MW49 Simulator Intravenous ARM III

M50B-B has been renovated to meet a wider range of training requirements.



Skills

- Anatomy of the arm
- Vacuum tube blood collection
- Palpation /Intravenous injection
- IV cannula placement and drip infusion
- Blood collection with a syringe/Flash-back confirmation

Set includes



1.	Arm model	qty: 2
2.	Supporting stand	qty: 2
3.	Injection pad (R pad)	qty: 2
4.	Injection pad (D pad)	qty: 2
5.	Injection pad	

(forearm/antebrachial) qty: 4 6. Injection pad(back of hand) qty: 2

7. Circulation pump with a plastic jar qty: 1

8. Bottle with stand (for sim blood) qty: 2

9. Simulated blood (10 swab pack) qty: 1 10. Arm vessels anatomy sleevelet qty: 2

11. Storage case qty: 1

12. Instruction manual qty: 1

Recommended devices


- IV cannula 23G ·syringe needle 23G
- ·butterfly needle 23G

Specifications

Size: W25×D11×H3cm each

Power: AC100-240V 50/60Hz Power Consumption: 5W

Four areas of puncture sites

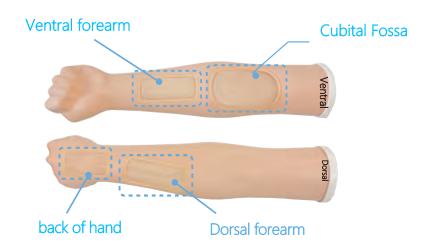

Two types of injection pads for venous line placement - works with bottles/stands





Pad for back of hand





Pad for ventral / dorsal of forearm



Two types of injection pads for cubital fossa - works with circulation pump





Regular pad

Realistic sensation. Texture, sensation and resistance to the needle. Three accessible venous lines represent basilic vein, median vein, and cephalic vein.

Durable pad

Durable training pad to maximize training opportunities. Long-life tissue

Both faces of the pad enclose two venous lines for training.

Consumables and replacement parts



R pad (set of 2) 11267-300-03



D pad (set of 2) 11430-020



pad for forearm/antebracheal (set of 2) pad for back of hand(set of 2) 11388-200




11388-300



supporting stand 11267-300-02



circulartion pump 11267-300-05

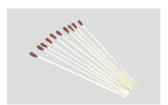

simulated blood (10 swab pack) 11388-400



arm vessels anatomy sleevelet 11430-010

KYOTO KAGAKU co..LTD



xw-kyoto@kyotokagaku.co.jp

Worldwide Inquires and Ordering

15 Kitanekoya-cho, Fushimi-ku,



Hong Kong Representative Office

Level 60. One Island East. 18 Westlands Road, Island East, Hong Kong





Head Office and Factories

Kyoto, 612-8388, Japan

European Representative Office

De-Saint-Exupery-Str.10, 60549 Frankfurt, Germany



+49-69-5060-28160 +49-69-5060-28170



All American regions KYOTOKAGAKU AMERICA INC.



info@kkamerica-inc.com

3109 Lomita Boulevard, Torrance, CA 90505-5108,USA



+1-310-325-8860



+1-310-325-8867